





# **HD Medical Video Recorder**

Record Full HD video from your endoscope or surgical camera. New state-of-the art recording platform records Full 1080p at 60 frames per second.



# Compact and robust

This rugged MVR features a stainless steel case and fluid resistant front panel. The compact size allows limitless mounting and installation options.



## **Intuitive Use**

It's so easy, no training is needed. Just plug it in and go. You'll be recording HD videos within minutes of opening the box.





# **Features**



#### Review

Review and control your photos and videos during the session or right before you close it.



### MVR Mount<sup>1</sup>

Flexible solutions for installation and cable management.



The included remote allows surgeons to record video and still images at the touch of a button.



#### **Patient Information**

Store patient data prior to recording for safe allocation. Easily use a Touch Screen Monitor<sup>2</sup> or a USB keyboard (optional).



#### Photos

Capture and store HD photos during your video recording.



#### **Foot Pedal**

Capture high quality HD video and HD photos to a USB drive using an endoscope trigger or foot pedal (optional).

### **Touch Screen Operation<sup>2</sup>**



Connect a USB Touch Screen Monitor and operate MVR Lite even more easily!

- <sup>1</sup> More information on our website
- <sup>2</sup> List of compatible Touch Screens on our website

# MediCap® MVR Lite Specs<sup>3</sup>

**HD Medical Video Recorder** Device type

> Records high-definition video and still images from medical video sources to

external USB drive.

Case/Buttons Stainless steel, fluid resistant

User interface Overlay on external monitor or touch

screen operation (record, storage

information)

**Video Inputs** HDMI, DVI (with adapter) Video Outputs HDMI, DVI (with adapter)

Video Formats H.264 / MPEG4

> Input video 1920x1200p 60 (reduced blanking only) 1080p 60/50/30/25 resolutions

1080i 60/50 1280x1024p50/60

720p 60/50

Image Format

Kit Contents

**Sound Recording** Mic by Stereo-Headset (4 poles) /

Line (not activated)

External Storage USB flash drive, USB hard drive

**USB Support** USB3.0, SuperSpeed, 1 front

connector, 2 rear connectors

Footswitch/ 3.5mm mini jack. Compatible with Trigger

standard footswitches (not included) and closed-contacts camera triggers

**Remote Control** IR remote control (included)

> Mechanical Size: 160 x 140 x 40mm /

6.3 x 5.5 x 1.6"

Weight: 0.76 kg / 1lb 11oz

**Environmental** Storage: -40° to +85° C

Operating: -20° to +35° C

Power AC 100-240V, 50/60Hz, 0.4-0.2A,

3-prong AC jack

**Certifications** CE 93/42/EEC, IEC 60601-1:2005,

ANSI/AAMI ES60601-1:2005, CAN/ CSA C22.2 No.60601, 1:2008, E361228,

FCC Part 15, FDA Class 1, RoHS

MediCap® MVR Lite. USB flash drive. HDMI cable, HDMI to DVI adapter, split adapter cable remote, power cable, remote control, Quick Start

Guide

<sup>3</sup> Subject to change



**Back Panel Connectors** 

#### MediCapture, Inc. -North America

2250 Hickory Road, Suite 200 Plymouth Meeting, PA 19462, USA



+1-888-922-7887


www.medicapture.com sales@medicapture.com

#### MediCapture Europe B.V.

Kraanpoort 10 6041 EG Roermond The Netherlands



(1) +31 475 214 305

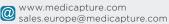

# MediCapture APAC Inc.

3FL-6, No. 432, Section 1, Keelung Road. Xinyi District 110, Taipei City, Taiwan, Republic Of China



+886-2-2345-3668



www.medicapture.com sales.asiapacific@medicapture.com